


InSight[®] Urinalysis

Veterinary Urine Analyser

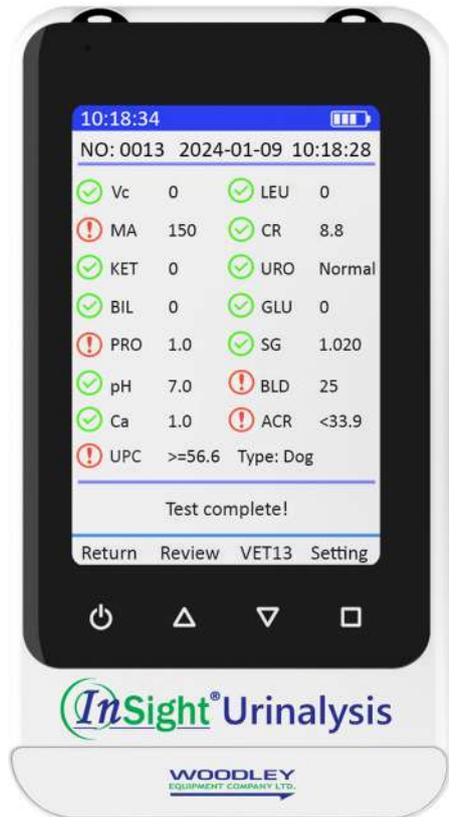

**13 PARAMETERS ON
ONE TEST STRIP
UPC RATIO ON EACH TEST**

"PASSIONATE ABOUT LABORATORY DIAGNOSTICS"

WOODLEY
EQUIPMENT COMPANY LTD.

InSight[®] Urinalysis

Veterinary Urine Analyser

Urinalysis is an important tool in disease detection, as well as monitoring and screening animal health. Abnormalities can be indicative of diseases of the urinary system as well as other organ systems, including liver function, acid base status and carbohydrate metabolism.

InSight Urinalysis is an easy to use Veterinary Urine Analyser that provides fast, accurate results and can test up to 13 parameters on one test strip.

FEATURES

- Easy to use
- Fast, accurate results
- Compact and portable
- Up to 13 parameters on one test strip
- Auto-recognition of test strips
- Stores up to 1,200 results
- Multiple species available
- Large LCD screen with touch keys
- Bluetooth connection to wireless printer
- Multi-language operation system
- Rechargeable battery
- Exclusive software for veterinary use

PARAMETERS

- Leucocytes
- Ketone
- Urobilinogen
- Bilirubin
- Glucose
- Protein
- Specific Gravity
- pH
- Blood
- Ascorbic Acid
- Microalbumin
- Calcium
- Creatinine

With Urine Protein/Creatinine Ratio
& Microalbumin/Creatinine Ratio

13 Parameters On
One Test Strip



SPECIFICATIONS

- Species:** Dog, Cat, Rabbit, Horse, Cattle, Sheep
- Measuring Principle:** Reflectance photometry
- Measurement Wavelength:** 470nm, 550nm, 620nm, 720nm
- Sample Supply Method:** Manual dipping
- Throughput:** Normal Mode - 70 tests per hour
Fast Mode - 140 tests per hour
- Measuring Mode:** Automatic single measuring mode
- Display:** LCD screen with touch keys
- Memory:** Up to 1,200 results
- Data Output:** USB, Bluetooth
- Operating Conditions:** 5°C ~ 40°C, ≤80% humidity
- Measuring Conditions:** 15°C ~ 30°C, ≤80% humidity
- Power Supply:** Rechargeable battery (300 tests/charge)
- Power Requirements:** 100V - 240V, 50/60Hz
- Dimensions:** 130 (L) x 70 (W) x 29 (H) mm
- Weight:** 200g



+44 (0) 1204 669033 (Option 1)



sales@woodleyequipment.com



www.woodleyequipment.com



Old Station Park Buildings
St. John Street
Horwich
Bolton
Lancashire
BL6 7NY
UK

**"PASSIONATE ABOUT
LABORATORY DIAGNOSTICS"**

**WOODLEY
EQUIPMENT COMPANY LTD.**