

PocketChem™ BA

Veterinary Blood Ammonia Measurement

World's only Handheld, Patient Side Blood Ammonia Analyser



IMMEDIATE TESTING

FAST, ACCURATE RESULTS.

IN THE PALM OF YOUR HAND.



"PASSIONATE ABOUT LABORATORY DIAGNOSTICS"

WOODLEY
EQUIPMENT COMPANY LTD.

PocketChem™ BA

Veterinary Blood Ammonia Measurement

Enables immediate testing and delivers fast, accurate results.
With its lightweight and compact design, blood ammonia measurement can be conducted in the palm of your hand

FEATURES

- Whole blood measurement
- Ready for use 24 hours a day
- Small 20µl sample volume
- Battery powered
- Lightweight and compact design
- Uses exclusive reagent strips without any pre-treatment
- Room temperature storage of reagent strips
- Detachable printer available



**World's only Handheld,
Patient Side
Blood Ammonia Analyser**

Available with detachable printer




APPLICATIONS/USES

- Supports the diagnosis of hepatic cephalopathy
- Helps detect and evaluate urea cycle defects
- Monitor patients with liver disease

SPECIFICATIONS

Sample	Whole blood
Measurement Item	Blood ammonia
Measurement Principle	Test strip, single wavelength reflection measurement method
Measurement Range	10-400N - µg/dL (8-285 µmol/L)
Wavelength	635nm (LED)
Sample Volume	20µl
Test Strip	Ammonia Test Kit II
Reaction Time	Approx. 3 minutes and 20 seconds/test
Warm-Up Time	1-2 minutes
Display	LCD
Operator Panel	Panel sheet (numeric and function keypad, 12 keys in total)

Memory Capacity	50 measurements
External Output	RS-232C
Temperature Compensation Function	Automatic compensation by built-in temperature sensor
Operating Temperature	10-35°C
Relative Humidity	25-85%
Power Supply	2 x AA-sized batteries or AC adapter
Dimensions	38 (H) X 125 (W) X 85 (D) mm 133mm (D with optional printer)
Weight	150g without batteries

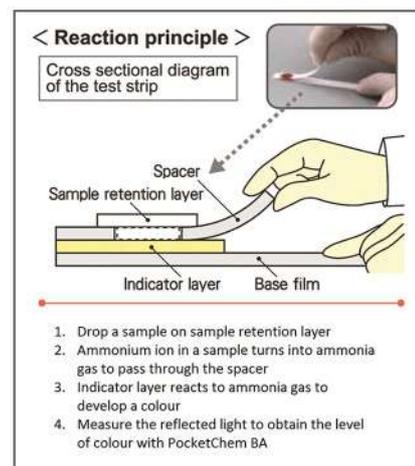

+44 (0) 1204 669033 (Option 1)

sales@woodleyequipment.com

www.woodleyequipment.com

Old Station Park Buildings
St. John Street
Horwich
Bolton
Lancashire
BL6 7NY
UK

**"PASSIONATE ABOUT
LABORATORY DIAGNOSTICS"**

**WOODLEY
EQUIPMENT COMPANY LTD.**